

14-5-2018



WESSELING BV

Quality since 1958

Rock 3 & 5

USER'S MANUAL AND MAINTENANCE GUIDE



CE

CONTENTS:

1. General information	2
- Manufacturer	- 2
- Equipment intended use and description	- 2
- storage	- 2
- Symbols	- 2
2. Special warnings and use restrictions	3
3. Use	3
- Starting up	- 3
- Setting up and adjustments	- 3
- Transportation and brakes	- 5
4. Cleaning and conservation	5
5. Maintenance	5
- Preventive maintenance	- 5
- Corrective maintenance	- 5
- Spare parts	- 5
6. Troubleshooting	6
7. Technical data sheet	7
8. Warranty	8
9. Service record	8

1. GENERAL INFORMATION

MANUFACTURER

Wesseling BV, Venenweg 31, 1161 AK Zwanenburg, Netherlands-Tel.+31 204975183-Fax+31 204975909 – website:www.wesseling-bv.com

EQUIPMENT INTENDED USE AND DESCRIPTION

The Rock tables are intended for examining and treating patients. With the additional options this table can be used for treatment and examination. The entire frame is epoxy powder coated and the top is upholstered and covered in fireproof synthetic leather. It's equipped with high quality actuators for height adjustment and electric flex position. Head section as well as back section can be inclined with the help of gas springs. The gas springs are all adjusted by handles that are easy to access on the side of the table or on the head section. The brakes and lock are centrally controlled by means of a foot pedal that is easy to activate.

STORAGE

When storing, the following conditions must be met:

1. Relative humidity 10% / 90 % - temperature -10°C / +50°C

SYMBOLS

The various symbols appearing in this manual and on the equipment itself are meant to draw the user's attention on important safety warnings and useful information. Please observe all warning, caution and hazard notices.

	WARNING Warning messages are used to advise users/patients about any risk of injury or danger if the specified procedures and conditions are not carefully observed during operation or maintenance activities.
	ELECTRIC SHOCK HAZARD This warning symbol identifies operating/maintenance instructions or procedures which, if not observed, could result in a risk of electric shock
	CAUTION Caution symbols are used to advise users about an operation or maintenance procedure, best practice or condition that may damage the equipment if not duly complied with.

2. SPECIAL WARNINGS AND USE RESTRICTIONS

	Equipment must not be used: • Without knowing the user's manual. • When the safety load exceeds the stated capacity. • Without the patient being monitored. • When the patient and/or equipment is not stable enough. • When table is not on its feet when stopped, during treatment, or when patient is moved to/from other equipment. • When the original product has been tampered with and/or modified from its original configuration. • When the controls and/or components are even only partially faulty or defective. • When incorrectly configured or when you think it may cause damage to a patient or user. • Do not allow any body parts to enter the areas with moving components and/or controls.
	• Do not connect to an electric system that does not comply with the current safety regulations. • Never spill liquid of any kind on the electric equipment. • Always unplug the mains power supply before working on any electrical parts or to move the equipment. • Do not damage the electric cables: avoid crushing; do not pull out the cable to disconnect, etc. • Electric cables should be checked regularly to see if they are damaged. In case of doubt, have them replaced.
	• Foam rubber and plastic components may burn if they come in contact with naked flames. • Corrosive and/or abrasive detergents may damage the covers and paintwork. • The equipment is not suitable for use outdoors and/or in surfaces that are not smooth enough. • UV and heat therapy lamps may damage the upholstery covering.
!!!!!!!!!!!! NEVER LEAVE THE TABLE UNATTENDED WHEN THERE IS A PATIENT ON IT !!!!!!!!!!!!!	

3. USE

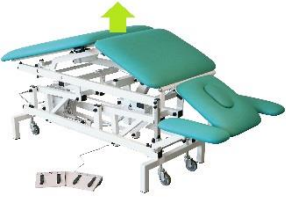
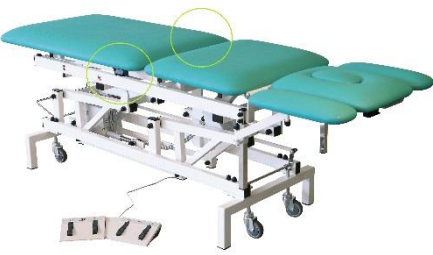
STARTING UP

For packaging reasons the couch is supplied covered by protective means that are to be removed, avoiding direct contact with cutters that could damage it. Any damage in transit must be reported to the transporter and seller immediately, without using the equipment for safety reasons.

SETTING UP AND ADJUSTMENTS

The table has an electric actuator for the height adjustment. The head section is fitted with a gas spring and the body section is mechanically adjustable.

ADJUSTMENT/FUNCTION	CONTROL	PROCEDURE
		Table is standard equipped with a 24v foot switch or hand switch.
		Push the relevant symbol to adjust the height of the entire table.

		To put the table in Pli-position(flex/drainage): Push the relevant symbol to reach the desired position WARNING • With the patient's weight on the section, loosen the lever only when you are sure that the entire load is supported. • Movement of the adjustable sections may cause pinching in the vicinity of the hinges.
		To adjust the head sections pull the control handle, and support the section as it moves. Release the handle to secure the position. WARNING • Do not sit on the head section • Movement of the adjustable sections may cause pinching in the vicinity of the hinges.
		To adjust the armrests loosen the lever of the armrest you want to adjust. Move the armrest in the desired direction and fasten the lever to keep it in place. WARNING • With the patient's weight on the section, loosen the lever only when you are sure that the entire load is supported. • Movement of the adjustable sections may cause pinching in the vicinity of the hinges. • Do not sit on the armrests.
		To adjust the sitting/leg section pull the control handle, and support the section as it moves. Release the handle to secure the position. The control handles are on both sides of the table. WARNING • With the patient's weight on the section, activate the control handle only when you are sure that the entire load is supported. • Movement of the adjustable sections may cause pinching in the vicinity of the hinges.

The table has four swivel wheels to allow movements in tight spaces. At the back-side of the table is the stirrup bar of the wheel-lift system. By pushing this stirrup bar downwards the table becomes mobile. Lift the stirrup bar upwards again after reaching the required position so the castors will be retracted. Never use the table if the castors are not fully retracted! The 4 feet can be adjusted individually when the floor is not even.

MOVING THE TABLE	PUT THE TABLE ON ITS FEET	PROCEDURE
		- To move the table use the foot to push down the bar till it stops. - To put the table on its feet use the foot to push up the bar until all the feet are on the floor.



WARNING



- Putting the table on its feet is essential when transferring and treating the patient.

4. USE AND SETTING UP OPTIONS

With the use of the optional accessories this table can be used for multiple treatments.

PAPH004: Paperrollholder Eco-Line



The paper holder is always mounted on the head-section of the Eco-Line. The maximum roll width that fits into holder 60cm.

Placing a roll of paper:

- Turn the knob and pull the rod out of the holder.
- Place the paper roll in the holder and push through the rod.
- Tighten the knob to hold the rod in place.

RONDOM3: Barswith system Eco-line



By choosing to barswith-system the standard switch will not be on the table. By pushing the bars which are located on both sides of the table, you can adjust the height of the table.

Adjusting the height:

- Push the bar left and right to adjust the height

Adjusting the flex position:

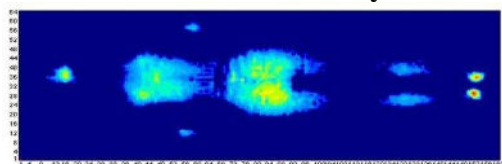
- Push the relevant symbol to reach the desired position

HSCH006+TSTUK00:
Extra handswitch


Standard the table is mounted with one switch of choice. We can equip the table with multiple switches. If you choose for an extra hand switch, you also need a link to connect two switches on the motor. If this option is chosen, the link (Y-cable) and the extra switch will be assembled on the table before delivery.

VSCH005+TSTUKROCK:
Extra footswitch

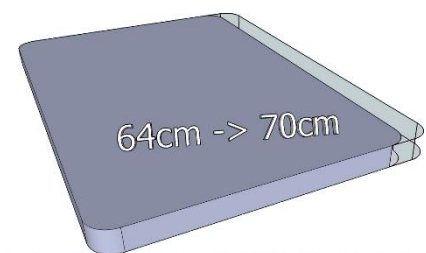

Standard the table is mounted with one switch of choice. We can equip the table with multiple switches. If you choose for an extra foot switch, you also need a link to connect two switches on the motor. If this option is chosen, the link (Y-cable) and the extra switch will be assembled on the table before delivery.

LAZYFOAM: memory foam


The standard parts are filled with polyether and are about 55mm thick. When the optional Lazy Foam is chosen the total thickness is 70mm. The upholstered sections will not only be thicker but also softer and therefore more comfortable for patients who need to be longer on the table.

MASK003: Neckcushion with velcro strap


This neck pillow ensures that the patient can sit/lie comfortably. The cushion has an elastic strap on the rear with Velcro closure. This is to position the cushion at almost every place. The cushion has a size of 58x18x9cm. The cover on the back is also equipped with a zipper and thus removable.

BREEDTE00: Different width 70cm


The standard width of a Praxis Elite is 64cm. If this is not wide enough, the table can be optionally expanded to 70cm. The length of the table will therefore not change.

DPUNIT00: Deskpanel under headsection

This panel is placed under the head section. In this way you can adjust the height of the table and the flexposition whilst sitting or standing on that side.

This part only works with TSTUCKROCK

5. CLEANING and CONSERVATION

For durable joy with the leather material, please clean regularly. For cleaning use a cloth or soft brush dampened using water and a neutral detergent. Any stains on the upholstery should be promptly removed with colourless ethyl alcohol diluted in 20% water. Rinse well.



WARNING



- ENSURE THE COUCH IS PERFECTLY CLEAN BY DISINFECTING THE SYNTHETIC LEATHER AND APPLYING A SINGLE-USE PAPER SHEET AT EVERY USE.
- HARSH DETERGENT PRODUCTS MAY DAMAGE THE SURFACES.
- DO NOT USE STEAM JETS, WASHING TUNNELS, WATER JET CLEANERS OR HIGH PRESSURE WASHERS

6. MAINTENANCE

Maintenance operations must be carried out by technicians who are familiar with the contents of this instruction booklet. Check the equipment is safe after any instances of "rough handling" (e.g. fall, bumps, etc.) and any time you are in doubt about safety.

ROUTINE MAINTENANCE

Routine maintenance operations must be carried out with the frequency indicated in the table.

FREQUENCY	ACTIVITY
Every 6 months or at least once a year	General equipment cleaning.
	Visual and functional inspection of the equipment and related hydraulic and manual controls.
	Frame and mechanical parts: Check screws and bolts are tight; state of wear of bushings and pins, presence of tears and deformations or damages.
	Oil leaks from the piston and/or gas springs must only be in the form of minor lubrication residue.
	Check that the adjustable sections are stable and do not give way when the control lever is released.
	Brake control: check regular functioning of the activating device and the efficiency of the brakes.
	Wheels: these must be securely fixed and in good condition, rolling smoothly and with their treads free of impediments.
	During use and operation, make sure there are no unusual sounds, vibrations or oscillations that may indicate actuator and/or mechanism malfunctions.
At the end of the activities above, complete the inspection performing a lift cycle with the stated maximum safe load.	

NON-ROUTINE MAINTENANCE

For safety reasons and to avoid invalidating the guarantee, no component may be opened or handled for repair. Any faults on the hydraulic piston, the gas springs and the wheel system must be diagnosed by Wesseling BV, who will organize the most suitable action.



WARNING



- IF SAFETY SEEMS TO BE AT RISK, DISABLE THE EQUIPMENT IMMEDIATELY PREVENTING ANY POSSIBILITY OF USE UNTIL SUCH RISK HAS BEEN REMOVED.
- PLEASE RECORD ANY DEFECTS/DAMAGES/SERVICES THAT ARE IMPORTANT FOR THE SAFETY OF THE EQUIPMENT IN THE APPROPRIATE FORM HEREWITH ENCLOSED.

SPARE PARTS

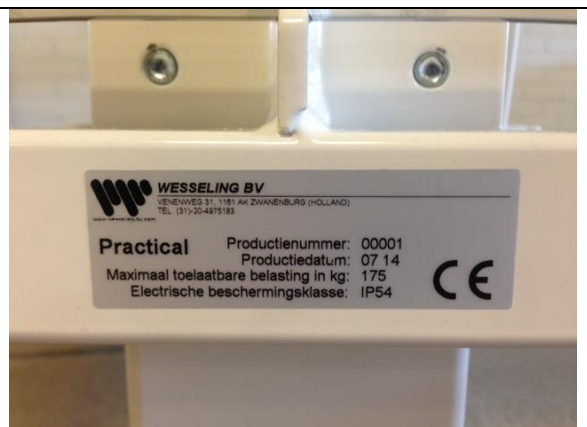
Only use original spare parts from Wesseling. Please contact Wesseling BV for the parts you might need.

7. TROUBLESHOOTING

DIRECT E-MAIL FOR TECHNICAL OR USERS ASSISTANCE:

info@wesseling-bv.com

If you rather want to speak to one of us on the phone you can call: +31-20 497 5183



Please note:

- Before contacting for service or help, please identify the "SERIAL NUMBER" in the product label attached to the equipment (see example above), or the transport document/invoice number. This will speed things up and improve service quality.

The equipment was tested in very stringent conditions. During testing no specific breakages or faults occurred.

If the equipment fails to operate, please follow the troubleshooting process below.

Issue	Cause	Solution
Adjustment of the head-/sitting section doesn't work or the position is not maintained	Gas spring faulty or incorrectly set	Contact Wesseling BV
	Transmission cable damaged or broken	
The actuator does not work	The attached switch is broken or not connected properly	Contact Wesseling BV
	Power cable is damaged or not connected properly	
	Actuator is damaged or not connected properly	

8. TECHNICAL DATA SHEET

Product name	Rock 3 & 5 Elektrisch
Product code	ROCK3E00 / ROCK5E00
Safe working load	200kg
Max patient weight	175kg
Height adjustment	54 tot 104 cm met een snelheid van 5,5mm/s - 6.000N, 470W, 2.0Ah, slag 200mm
Adjustment system	Actuator uit Denemarken
Adjusting electrics	24v voetschakelaar
Sections adjusted by gas spring	2
Sections adjusted electric	1
Head section	-60° to +75°
Back section	0° to +90°
Middle section	0° to +45°
Adjustment Headsection	By lever under the headsection
Adjustment back section	By levers on both sides of the table
Adjustment middle section	By footswitch
Casters	Ø75mm
Braking system	Central braking system at the back
Equipment weight	± 90kg
Length	200cm
Width	64cm or 70cm (optional)
Upholstery thickness	± 55mm
Filling	Polyether SG 90/30
Semi Leather	Standard artificial leather, Original SKAI® or STAMSKIN
Dimensions	