



## Praxis Elite

User Manual and Maintenance Instructions



Doc. Number: IFU- EN - PRAXIS-001

Version: 1.0

Date: 17-03-2026



# Table of Contents

|           |                                                            |                                            |
|-----------|------------------------------------------------------------|--------------------------------------------|
| <b>1.</b> | <b>GENERAL INFORMATION.....</b>                            | <b>4</b>                                   |
| 1.1       | <i>Purpose of this document.....</i>                       | 4                                          |
| 1.2       | <i>Manufacturer .....</i>                                  | 4                                          |
| 1.3       | <i>Product status and regulations.....</i>                 | 4                                          |
| 1.4       | <i>User profile .....</i>                                  | 4                                          |
| 1.5       | <i>Symbols.....</i>                                        | 4                                          |
| <b>2.</b> | <b>Application and context of use.....</b>                 | <b>5</b>                                   |
| 2.1       | <i>Intended use .....</i>                                  | 5                                          |
| 2.2       | <i>Users and patients.....</i>                             | 5                                          |
| 2.3       | <i>Application environment .....</i>                       | 5                                          |
| 2.4       | <i>Use limitations.....</i>                                | 5                                          |
| 2.5       | <i>Storage and transport .....</i>                         | 5                                          |
| 2.6       | <i>Manufacturer's declarations .....</i>                   | 5                                          |
| <b>3.</b> | <b>Safety and warnings.....</b>                            | <b>6</b>                                   |
| 3.1       | <i>General safety principles .....</i>                     | 6                                          |
| 3.2       | <i>Use-related risks .....</i>                             | 6                                          |
| 3.3       | <i>Safe use .....</i>                                      | 6                                          |
| 3.4       | <i>Incidents and reporting .....</i>                       | 7                                          |
| <b>4.</b> | <b>Product design and configuration.....</b>               | <b>8</b>                                   |
| 4.1       | <i>General configuration.....</i>                          | 8                                          |
| <b>5.</b> | <b>Technical specifications.....</b>                       | <b>13</b>                                  |
| 5.1       | <i>Technical data .....</i>                                | 13                                         |
| <b>6.</b> | <b>Installation and commissioning.....</b>                 | <b>14</b>                                  |
| 6.1       | <i>Delivery.....</i>                                       | 14                                         |
| 6.2       | <i>Inspection upon receipt .....</i>                       | 14                                         |
| 6.3       | <i>First use .....</i>                                     | 14                                         |
| <b>7.</b> | <b>Use and operation.....</b>                              | <b>15</b>                                  |
| 7.1       | <i>Pre-use conditions .....</i>                            | 15                                         |
| 7.2       | <i>Operation.....</i>                                      | 15                                         |
| 7.2.1     | <i>Height adjustment .....</i>                             | 15                                         |
| 7.2.1.1   | <i>Electric .....</i>                                      | 15                                         |
| 7.2.1.2   | <i>Hydraulic .....</i>                                     | 15                                         |
| 7.2.2     | <i>Adjusting the upholstery sections .....</i>             | <b>Fout! Bladwijzer niet gedefinieerd.</b> |
| 7.2.2.1   | <i>Gas spring operation.....</i>                           | 16                                         |
| 7.2.2.2   | <i>Electrically adjustable .....</i>                       | 16                                         |
| 7.2.2.3   | <i>Face opening.....</i>                                   | 16                                         |
| 7.2.3     | <i>Moving and positioning .....</i>                        | 16                                         |
| 7.3       | <i>Accessories &amp; additional options.....</i>           | 18                                         |
| 7.3.1     | <i>Available accessories &amp; additional options.....</i> | 18                                         |
| 7.3.2     | <i>Using accessories .....</i>                             | 19                                         |

---

|            |                                         |           |
|------------|-----------------------------------------|-----------|
| 7.3.2.1    | Paper roll holders.....                 | 19        |
| 7.3.2.2    | Adapters.....                           | 19        |
| 7.3.2.3    | Supports.....                           | 20        |
| 7.3.3      | Using additional options.....           | 22        |
| 7.3.3.1    | Side rails.....                         | 22        |
| 7.3.3.2    | Raising kit.....                        | 23        |
| 7.3.3.3    | Trapeze bar.....                        | 23        |
| <b>8.</b>  | <b>Cleaning and hygiene .....</b>       | <b>24</b> |
| 8.1        | <i>Cleaning</i> .....                   | 24        |
| 8.2        | <i>Disinfection</i> .....               | 24        |
| 8.3        | <i>Restrictions</i> .....               | 24        |
| <b>9.</b>  | <b>Maintenance and inspection .....</b> | <b>25</b> |
| 9.1        | <i>Daily inspection</i> .....           | 25        |
| 9.2        | <i>Preventive maintenance</i> .....     | 25        |
| <b>10.</b> | <b>Faults and support.....</b>          | <b>26</b> |
| 10.1       | <i>General</i> .....                    | 26        |
| 10.2       | <i>Common faults</i> .....              | 26        |
| 10.3       | <i>Support</i> .....                    | 26        |
| <b>11.</b> | <b>Disposal.....</b>                    | <b>27</b> |
| 11.1       | <i>General</i> .....                    | 27        |
| 11.2       | <i>Reuse and circularity</i> .....      | 27        |
| <b>12.</b> | <b>Warranty .....</b>                   | <b>27</b> |
| 12.1       | <i>Warranty terms</i> .....             | 27        |
| 12.2       | <i>Exclusions</i> .....                 | 27        |
| 12.3       | <i>Repairs</i> .....                    | 27        |
| 12.4       | <i>Spare parts</i> .....                | 28        |

---

## 1. GENERAL INFORMATION

### 1.1 Purpose of this document

This manual provides information for the safe, correct and practical use of the device and is intended for both specialist personnel and end users. We recommend reading the entire manual carefully before putting the product into use.

For questions or clarification, please contact your supplier, who can provide appropriate advice.

The importance of consulting and understanding this manual is indicated on the product by the following symbol:



**Follow the instructions for use**

### 1.2 Manufacturer

**Wesseling BV**  
**Venenweg 31**  
**1161 AK Zwanenburg,**  
**The Netherlands**  
**Tel. +31 20 497 5183**  
**Website: [www.wesseling-bv.com](http://www.wesseling-bv.com)**  
**Email: [info@wesseling-bv.com](mailto:info@wesseling-bv.com)**

### 1.3 Product status and regulations

This product is a Class I medical device in accordance with Regulation (EU) 2017/745 (MDR). The product complies with the general safety and performance requirements set out in Annex I.

### 1.4 User profile

The product is intended exclusively for professional use by qualified personnel.

### 1.5 Symbols

The symbols used in these instructions and on the equipment are intended to draw the user's attention to important safety warnings and useful information. Observe all warnings and restrictions and remain alert to all indicated hazards.

#### Note!

“Note!” indicates recommended actions concerning use, maintenance and the condition of the table where there is no immediate risk of damage or injury to the user. Failure to observe these notes may nevertheless damage the equipment.




#### Warnings

Warning notices alert users and patients to a risk of injury or danger if the specified procedures and conditions are not carefully observed during operation or maintenance activities.



### **Danger**

This warning symbol identifies maintenance instructions or procedures that, if not observed, may result in a risk of electric shock.
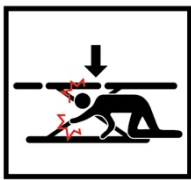
|                                                                                                                                                                                               |                                                                                                                                                                                    |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|                             |               |
| Warning: risk of entrapment or injury when adjusting the treatment table. Keep hands, head and other body parts away from moving parts and from beneath the patient surface during operation. | Warning: risk of entrapment when adjusting the treatment table electrically or hydraulically. Keep hands, feet and legs away from the scissor mechanism and base during operation. |

## **2. Application and context of use**

### **2.1 Intended use**

The treatment table is intended to support and position patients or clients during medical, therapeutic or cosmetic treatments.

The product enables a stable and ergonomic working posture.

### **2.2 Users and patients**

Use is permitted only by professionals. Patients must not operate the product independently.

### **2.3 Application environment**

Suitable for use in professional indoor environments such as treatment rooms and medical practices.

### **2.4 Use limitations**

Use is not permitted:

- Outside professional indoor environments
- Beyond the specified load limits
- For purposes other than those described

### **2.5 Storage and transport**

- Temperature: -10 °C to +50 °C
- Relative humidity: 10% to 90%

### **2.6 Manufacturer's declarations**

- The product is not a measuring instrument
- The product is not intended for clinical investigation
- The product is not sterile and is not intended for single use
- Correct operation requires maintenance in accordance with these instructions for use

- Use is limited to the applications described

## 3. Safety and warnings

### 3.1 General safety principles

Before each use, check the product for correct and safe operation.

### 3.2 Use-related risks

| Risk code | Warning                                                                                |
|-----------|----------------------------------------------------------------------------------------|
| R-01      | <b>Entrapment hazard: Keep hands, feet and other body parts away from moving parts</b> |
| R-02      | <b>Electrical hazard: Do not use if cables are damaged</b>                             |
| R-03      | <b>Tipping hazard: Ensure a stable surface</b>                                         |

#### Note!

*The following residual risks remain despite the measures taken:*

- *Entrapment at moving parts*
- *Electrical hazard due to defective components*

### 3.3 Safe use

To use the treatment table safely, the user must strictly observe the instructions below. Failure to follow these instructions may result in injury to the patient or user, or damage to the product.

The treatment table must not be used in the following situations:

- If the user has not read and fully understood the instructions for use
- If the load exceeds the specified maximum capacity
- Without continuous visual contact with the patient
- If the patient and/or table is not positioned securely
- If the table is not resting fully on its feet while stationary, during treatment or during patient transfer
- If the product has been adapted or modified from its original configuration
- If controls or components are defective or partially defective
- If settings are incorrect or could pose a risk to the patient or user



#### Warnings

*Use the treatment table only when it is correctly adjusted and stable, to prevent falls and injury.*

The following safety instructions must also be observed during use:

- Keep body parts away from moving components and hinge points
- Prevent entrapment when adjusting components
- Ensure that the patient remains correctly positioned during every movement



#### Warnings

---

*Incorrect use of adjustable sections may cause entrapment or sudden movement of the patient.*

The following additional instructions apply to products with electrical components:

- Connect the table only to an electrical installation that complies with the applicable safety standards
- Prevent liquids from coming into contact with electrical components
- Always disconnect the power supply before performing maintenance or moving the table
- Prevent damage to cables (do not pull, kink or trap them)
- Regularly inspect cables for wear or damage and replace them if necessary



**Danger**

*Never work on electrical components while the product is connected to the power supply.*

### **3.4 Incidents and reporting**

Serious incidents must be reported to the manufacturer.

## 4. Product design and configuration

### 4.1 General configuration

The product is available in various configurations that differ in:

- Frame colour
- Height adjustment and controls
- Versions
- Width and openings

The product code shown on the product label identifies all features of the device configuration.

Use the code below to quickly identify your product configuration.

|   |   |   |   |   |   |   |   |   |   |    |    |   |    |    |    |
|---|---|---|---|---|---|---|---|---|---|----|----|---|----|----|----|
| 0 | 1 | * | * | . | * | * | * | * | * | *  | *  | . | *  | *  | *  |
| 1 | 2 | 3 | 4 | . | 5 | 6 | 7 | 8 | 9 | 10 | 11 | . | 12 | 13 | 14 |

#### Position 1-2

Indicates the model name: 01 = Praxis

#### Position 3-4

Indicates the model type

|                    |    |
|--------------------|----|
| (Praxis) 1 Junior  | 01 |
| (Praxis) 1         | 02 |
| (Praxis) 2         | 03 |
| (Praxis) 3         | 04 |
| (Praxis) 3 Comfort | 05 |
| (Praxis) 4         | 06 |
| (Praxis) 4 Comfort | 07 |
| (Praxis) 3 Podo    | 08 |
| (Praxis) 4 Podo    | 09 |

|                                                                                     |                                                                                     |                                                                                      |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
|  |  |  |
| 01 + 02                                                                             | 03                                                                                  | 04                                                                                   |
|  |  |  |
| 05                                                                                  | 06                                                                                  | 07                                                                                   |

|                                                                                   |                                                                                    |
|-----------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
|  |  |
| 08                                                                                | 09                                                                                 |

### Position 5

Indicates the frame colour:

|                 |   |
|-----------------|---|
| White, RAL 9010 | 1 |
| Black, RAL 9005 | 2 |

|                                                                                   |                                                                                   |
|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|
|  |  |
| 1                                                                                 | 2                                                                                 |

### Note!

*Only a white frame is available for the Praxis range.*

### Position 6

Indicates the height adjustment type:

|                             |   |
|-----------------------------|---|
| Electric height adjustment  | 1 |
| Hydraulic height adjustment | 2 |

|                                                                                     |                                                                                      |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
|  |  |
| 1                                                                                   | 2                                                                                    |

## Position 7+8

Indicates the control type:

|                                 |    |
|---------------------------------|----|
| Pump bar                        | 00 |
| Foot switch                     | 01 |
| Foot switch in tray holder      | 02 |
| All-round switch                | 03 |
| Hand control                    | 04 |
| Dual hand control               | 05 |
| Dual foot switch in tray holder | 06 |
| Dual foot switch                | 07 |

|                                                                                    |                                                                                    |                                                                                     |                                                                                      |
|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
|   |   |   |   |
| 00                                                                                 | 01                                                                                 | 02                                                                                  | 03                                                                                   |
|  |  |  |  |
| 04                                                                                 | 05                                                                                 | 06                                                                                  | 07                                                                                   |

## Position 9

Control box with Safety Anti-Collision software (SAC System).

Function:

|     |   |
|-----|---|
| No  | 0 |
| Yes | 1 |

|                                                                                     |                                                                                      |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
|  |  |
| 0                                                                                   | 1                                                                                    |

### Position 10

Control box with 250 kg maximum lifting capacity (MAXCAP250)

Function: Increases maximum lifting capacity from 200 kg to 250 kg.

|     |   |
|-----|---|
| No  | 0 |
| Yes | 1 |



0



1

### Position 11

Second actuator for electric backrest adjustment.

#### Note!

NOT available on: Praxis 1 JR/Praxis 1/Praxis 3 Comfort/Praxis 4 Comfort

|     |   |
|-----|---|
| No  | 0 |
| Yes | 1 |



0



1

### Position 12

Upholstery width

|       |   |
|-------|---|
| 64 cm | 1 |
| 70 cm | 2 |
| 80 cm | 3 |

### Position 13

Face opening

|     |   |
|-----|---|
| No  | 0 |
| Yes | 1 |



0



1

---

**Position 14**

Opening in the seat section.

**Note!**

*Applicable to Kombi (S) only.*

|     |   |
|-----|---|
| No  | 0 |
| Yes | 1 |

|                                                                                   |                                                                                     |
|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|
|  |  |
| 0                                                                                 | 1                                                                                   |

## 5. Technical specifications

### 5.1 Technical data

|                              | Praxis 1 JR                                                                 | Praxis 1 | Praxis 2    | Praxis 3     | Praxis 4     |
|------------------------------|-----------------------------------------------------------------------------|----------|-------------|--------------|--------------|
| Safe working load            | 400 kg                                                                      |          |             |              |              |
| Maximum patient weight       | 200   250 kg                                                                |          |             |              |              |
| Electric height adjustment   | 42 to 112 cm at a speed of 7.5 mm/s - 8,000 N, 470 W, 2.0 Ah, stroke 300 mm |          |             |              |              |
| Hydraulic height adjustment  | 48 to 98 cm, hydraulic pump, 10,000 N, stroke 200 mm                        |          |             |              |              |
| Adjustable sections          | -                                                                           | -        | 1           | 2            | 3            |
| Adjustable head/back section | -                                                                           | -        | 0° tot +85° | 0° tot +85°  | 0° tot +85°  |
| Adjustable foot section      | -                                                                           | -        |             | -85° to +35° | -85° to +35° |
| Castors                      | Ø75mm                                                                       |          |             |              |              |
| Brake system                 | Central control at the rear of the table                                    |          |             |              |              |
| Weight                       | ±100kg                                                                      | ±100kg   | ±100kg      | ±100kg       | ±100kg       |
| Length                       | 198 cm                                                                      | 198 cm   | 198 cm      | 198 cm       | 202 cm       |
| Width                        | 64 cm   70 cm   80 cm                                                       |          |             |              |              |
| Upholstery thickness         | ± 60 mm   ± 70 mm                                                           |          |             |              |              |
| Padding                      | Polyether SG 90/30                                                          |          |             |              |              |
| Artificial leather           | Fresh, Purfect, Tovega                                                      |          |             |              |              |

|                              | Praxis 3 Comfort                                                            | Praxis 4 Comfort | Praxis 3 Podo                                                               | Praxis 4 Podo |
|------------------------------|-----------------------------------------------------------------------------|------------------|-----------------------------------------------------------------------------|---------------|
| Safe working load            | 400 kg                                                                      |                  |                                                                             |               |
| Maximum patient weight       | 200   250 kg                                                                |                  |                                                                             |               |
| Electric height adjustment   | 60 to 130 cm at a speed of 7.5 mm/s - 8,000 N, 470 W, 2.0 Ah, stroke 300 mm |                  | 42 to 112 cm at a speed of 7.5 mm/s - 8,000 N, 470 W, 2.0 Ah, stroke 300 mm |               |
| Hydraulic height adjustment  | 66 to 116 cm, hydraulic pump, 10,000 N, stroke 200 mm                       |                  | 48 to 98 cm, hydraulic pump, 10,000 N, stroke 200 mm                        |               |
| Adjustable sections          | 2                                                                           | 3                | 3                                                                           | 4             |
| Adjustable head/back section | 0° tot +85°                                                                 |                  |                                                                             |               |
| Adjustable foot section      | -90° to +35°                                                                |                  | -85° to +35°                                                                |               |
| Castors                      | ø75mm                                                                       |                  |                                                                             |               |
| Brake system                 | Central control at the rear of the table                                    |                  |                                                                             |               |
| Weight                       | ±100kg                                                                      | ±100kg           | ±100kg                                                                      | ±100kg        |
| Length                       | 198 cm                                                                      | 202 cm           | 198 cm                                                                      | 202 cm        |
| Width                        | 64 cm   70 cm   80 cm                                                       |                  |                                                                             |               |
| Upholstery thickness         | ± 60 mm   ± 70 mm                                                           |                  |                                                                             |               |
| Padding                      | Polyether SG 90/30                                                          |                  |                                                                             |               |
| Artificial leather           | Fresh, Purfect, Tovega                                                      |                  |                                                                             |               |

---

## 6. Installation and commissioning

### 6.1 Delivery

This product is delivered fully assembled and is ready for immediate use. During transport, the upholstery is protected by cardboard packaging that can be removed easily.

Export shipments are additionally packed in a suitable box for protection during transport.

#### Note!

*When removing the packaging, ensure that scissors and/or knives do not come into direct contact with the product. They may damage the upholstery!*

### 6.2 Inspection upon receipt

Any damage sustained during transport must be reported immediately to the carrier and supplier. For safety reasons, do not use the equipment.

### 6.3 First use

Check all functions before first use.

Hydraulic versions may be used immediately after inspection. For electric versions, first release the power cord and control secured with a cable tie. The power cord can then be plugged into the mains socket.



#### Warnings



#### Danger

*See Chapter 7, "Use and operation"*

## 7. Use and operation

### 7.1 Pre-use conditions

Carry out the following checks before each use:

1. Inspect the table for visible damage to items such as the upholstery, frame and cables
2. Check the stability of the table. On an uneven floor, adjust the levelling feet beneath the table legs. See illustration.
3. Check that the power cable is intact and correctly connected



*Turn the foot cap to adjust it to the required height.*



### Warnings

*Do not use the table if defects or irregularities are found*

### 7.2 Operation

#### 7.2.1 Height adjustment

##### 7.2.1.1 Electric

Operate the table using the foot, hand or perimeter control.



Foot control



Hand control



Perimeter control

##### 7.2.1.2 Hydraulic

Operate the table using the pump bar.



Hydraulic



## Warnings

Keep hands, feet and cables away from moving parts. (See R-01)

### 7.2.2 Adjusting the upholstery sections


#### 7.2.2.1 Gas spring operation

Use the gas spring levers and support the section while adjusting it.



#### 7.2.2.2 Electrically adjustable

Use a hand or foot control to adjust the back section electrically.

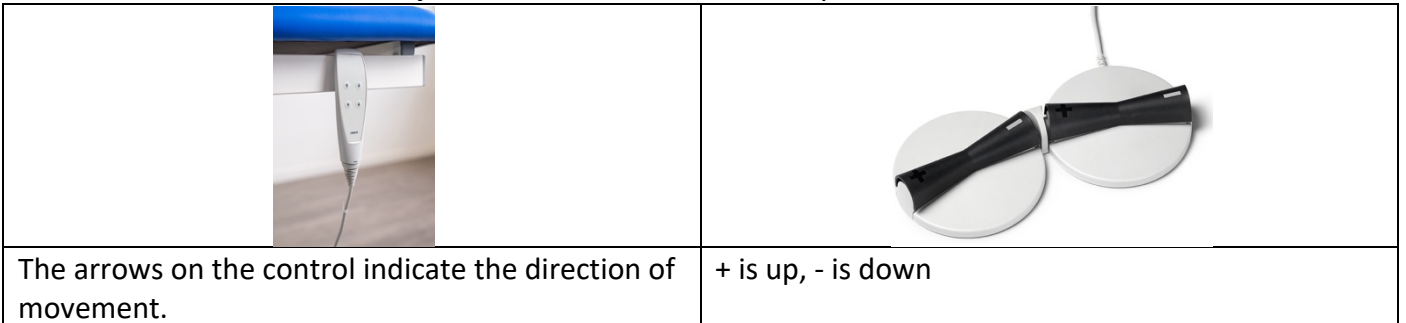

#### 7.2.2.3 Face opening

Almost every table can be supplied with a face opening. When this option is selected, a matching insert in the table colour is supplied as standard.

#### Note!

Not available for the Praxis 1 JR/Praxis 1



### 7.2.3 Moving and positioning

Praxis tables have four twin-wheel castors for manoeuvring in confined spaces. The wheel-lift system bar is located at the rear of the table. Pressing this bar down makes the table mobile. Lift the bar after reaching the required position so that the castors retract again. Never use the table unless the castors are fully retracted! The four legs have levelling feet that can be adjusted individually when the floor is uneven.

#### Note!

For adjusting the feet, see 7.1 Pre-use conditions.

| MOVING THE TABLE                                                                  | LOWERING THE TABLE ONTO ITS FEET                                                  | PROCEDURE                                                                                                                                                                                                                                     |
|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  |  | <ul style="list-style-type: none"> <li>- To move the table, press the bar down with your foot until it stops.</li> <li>- To lower the table onto its feet, push the bar up with your foot until all feet are resting on the floor.</li> </ul> |



### Warnings

*The table must be resting on its feet during patient transfer or treatment.*

## 7.3 Accessories & additional options

### 7.3.1 Available accessories & additional options

| Item code                  | Praxis 1<br>JR | Praxis<br>1 | Praxis<br>2 | Praxis<br>3 | Praxis<br>4 | Praxis 3<br>Comfort | Praxis 4<br>Comfort | Praxis<br>3 Podo | Praxis<br>4 Podo |
|----------------------------|----------------|-------------|-------------|-------------|-------------|---------------------|---------------------|------------------|------------------|
| VERHOG00                   | X              | X           | X           | X           | X           | X                   | X                   | X                | X                |
| HEKKEN00                   | X              | X           | X           | -           | -           | -                   | -                   | -                | -                |
| HEKKEN01                   | -              | X           | -           | -           | -           | -                   | -                   | -                | -                |
| PAPE003                    | -              | -           | X           | -           | -           | -                   | -                   | -                | -                |
| PAPRAX                     | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| ADAPB01                    | -              | -           | X           | -           | -           | -                   | -                   | -                | -                |
| ADAPB02                    | -              | -           | -           | X           | X           | -                   | -                   | X                | x                |
| ADAPB03From<br>80 cm width | -              | -           | -           | X           | X           | -                   | -                   | X                | X                |
| VOET001                    | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| HAKST00                    | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| BEENS01                    | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| BEENS02                    | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| ARMST00                    | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| ARMSTSAAR3D                | -              | -           | X           | X           | X           | X                   | X                   | X                | X                |
| ARMST3                     | -              | -           | -           | X           | X           | -                   | -                   | X                | X                |

#### Note!

*Supports can be fitted only in combination with adapters!*

### 7.3.2 Using accessories

#### 7.3.2.1 Paper roll holders

Item code: PAPRAX



On Praxis tables, the paper roll holder is always mounted on the head/back section.

Fitting a paper roll:

1. Remove the stainless-steel rod from the holder.
2. Push the rod through the roll.
3. Place it back in the holder.

#### Note!

*The maximum roll width that fits in the holder is 60 cm.*

#### 7.3.2.2 Adapters

The adapters are supplied as a set and fitted to the table before delivery. The adapters are adjustable, allowing the various supports to be positioned correctly.

Adjust the height of the relevant support as follows:

1. Loosen the black clamping handle and move the support to the required height.
2. Then tighten the clamping handle securely to keep the support at the selected height.



#### Warnings

*When the patient's weight is on the support, loosen the clamping handle only after ensuring that the full weight is supported!*

Item code: ADAPB02 + ADAPB03 + ADAPB05 + ADAPB06

These adapters can be locked in any position using the pull knob.



Item code: ADAPB01 + ADAPB04



### 7.3.2.3 Supports

#### Foot supports

The foot supports are supplied as a set and delivered separately with the table. They can be positioned inward or outward. In combination with the swivel adapters, the supports can always be positioned correctly.

Item code: VOET001



#### Heel supports

The heel supports are supplied as a set and delivered separately with the table. In combination with the swivel adapters, they can always be positioned correctly.

Item code: HAKST00



### Leg supports

The leg/knee supports are supplied as a set and delivered separately with the table. The angled upright, in combination with the adapters, allows the supports to be adjusted in width, depth and height, ensuring correct positioning.

**Item code: BEENS01**



### Arm support

The arm supports are supplied individually and delivered separately with the table. The angled upright, in combination with the adapters, allows the arm supports to be adjusted in width, depth and height, ensuring correct positioning.

**Item code: ARMST00**



### 3D arm support

This PU arm support is supplied individually and can be adjusted to any angle. Turn the ball knob on the underside to set the required angle. The angled upright, in combination with the adapters, allows the support to be adjusted in width, depth and height, ensuring correct positioning. The arm pad measures 350 x 115 mm and is ergonomically shaped.

**Item code: ARMSTSAAR3D**



### Warnings

*When the patient's weight is on the support, loosen the clamping handle only after ensuring that the full weight is supported!*

### Folding arm support

These arm supports are permanently mounted on the table and can be folded backwards. They are supplied as standard in the same colour as the table.

**Item code: ARMST3**



#### Note!

Available *ONLY* for the following models: Praxis 3/Praxis 4/Praxis 3 Podo/Praxis 4 Podo



#### Warnings

*Do not use the supports as push-off handles when standing up. The support could tip backwards.*

### 7.3.3 Using additional options

#### 7.3.3.1 Side rails

The optional side rails provide additional support and give the patient greater confidence when getting on and off the table and during treatment. The side rails can be folded up and down easily and must be checked for correct locking before each use.



#### Warnings

*Never leave a patient unattended, even when the side rails are raised. The side rails provide additional support but do not guarantee protection against falls or other unsafe situations. Always keep the patient under supervision.*

**Item code:**

**HEKKEN00:** 100 cm long

**HEKKEN02:** 140 cm long (not available for Praxis)



### 7.3.3.2 Raising kit

The optional raising kit raises the treatment table's lower frame by 6 cm. This increases the total clearance to 17 cm, making it easier to position and move aids such as patient lifts beneath the treatment table.

#### Note!

*When combined with a hydraulic pump, the pump operating lever is also positioned higher.*




#### Warnings

*Always check that the treatment table is stable before use.*

**Item code: VERHOG00**

This raises the frame by 6 cm.

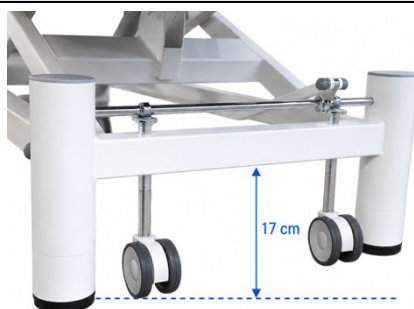

### 7.3.3.3 Trapeze bar

The trapeze bar is attached to the upper frame, so the distance between the trapeze bar and the patient remains constant when the height is adjusted. It is supplied as standard with a height-adjustable triangular handgrip.

#### Note!

- *Available only for the Praxis 2.*
- *The total length of the table will increase by 7.5 cm.*

**Item code: PAPE003**



---

## 8. Cleaning and hygiene

Proper cleaning and disinfection are essential for preserving the upholstery and maintaining a hygienic working environment. Follow the instructions below carefully.

### 8.1 Cleaning

Clean the table regularly to extend the service life of the upholstery.

#### PVC artificial leather (Skai FRESH)

Clean with a soft cloth or brush, water and a neutral cleaning agent. Remove stains immediately. Always rinse with clean water and dry the surface thoroughly.

#### PU artificial leather (Skai PERFECT)

Clean with a soft cloth or brush and a cleaning agent diluted with 70% water. Always rinse with clean water and dry the surface thoroughly.

#### Note!

*A list of tested and approved cleaning agents and disinfectants for each upholstery type is available on request. Colour charts with the corresponding specifications are also available on request.*

### 8.2 Disinfection

Ensure that the table is completely clean after each use. Disinfect the artificial leather and, where appropriate, use a new paper roll or protective cover for each new patient.



#### Warnings

*Use only suitable disinfectants that are compatible with the artificial leather. Excessive or incorrect use may shorten the service life of the upholstery.*

### 8.3 Restrictions

- Do not use aggressive cleaning agents such as solvents, chlorine or abrasive products
- Do not use steam cleaners or pressure washers.
- Avoid prolonged exposure to strong chemicals.
- Remove stains as quickly as possible to prevent permanent damage.



#### Warnings

*Failure to follow these guidelines may damage the upholstery and void the warranty!*

## 9. Maintenance and inspection

### 9.1 Daily inspection

Before each use, the user must check that the product operates safely and correctly.

#### Note!

*Check that the functions operate as described in 7.1 Pre-use conditions.*



#### Warnings

*Do not use the product if:*

- *Components are loose or damaged.*
- *The table is not stable.*
- *Safe treatment cannot be guaranteed*

### 9.2 Preventive maintenance

Periodic maintenance is required for safe, long-term operation.

#### Recommended frequency:

- At least once a year, depending on intensity of use

#### Preventive maintenance may be carried out by:

- The user, for visual and functional checks.
- The manufacturer or a party authorised by the manufacturer, for comprehensive inspection and maintenance.

#### Components to be checked:

- Frame structure for cracks or deformation
- Fasteners such as bolts and screws
- Operation of gas springs and adjustment mechanisms
- Wheel-lift system and castor mobility
- Electrical components and wiring (if applicable)

#### Note!

*During use, watch for:*

- *Unusual noises*
- *Jerky movements*
- *Reduced stability*

*These may indicate wear or defects.*

- *Failure to perform maintenance correctly may reduce the safety and operation of the product.*
- *If there is any doubt about the technical condition of the product, maintenance must be performed by the manufacturer or an authorised party.*



#### Warnings

- *Do not perform maintenance while a patient is on the table.*
- *Ensure that the product is stable and cannot move during maintenance work.*



#### Danger

*Disconnect the power supply before working on electrical components (if applicable)*

## 10. Faults and support

### 10.1 General

In the event of a fault or abnormal operation:

- Stop using the product
- Check whether the cause is visible (e.g. cable, control, obstruction)
- Do not use the product if safety cannot be guaranteed.



#### Warnings

*Stop using the product if a fault is identified that could affect the safety of the patient or user.*

### 10.2 Common faults

**Problem: Height adjustment does not work**

**Possible causes:**

- Control is defective/not connected correctly
- Power cord is defective/not connected correctly
- Actuator is defective/not connected correctly
- Control box is defective/not connected correctly
- Hydraulic pump bars have become detached from the pump

**Problem: An adjustable section cannot be positioned or will not remain in position.**

**Possible causes:**

- Defective or incorrectly adjusted gas spring
- Cable is damaged or broken



#### Warnings

*Do not use the product if adjustable sections move unexpectedly or fail to remain in position, to prevent entrapment or injury.*

#### Note!

*Check simple causes such as the power supply and controls before contacting the manufacturer.*

### 10.3 Support

For technical support, contact the manufacturer. See 1.2 Manufacturer.

#### Note!

*Have the product serial number available when contacting support.*



## 11. Disposal

### 11.1 General

At the end of its service life, the product must be disposed of in accordance with applicable local regulations.

#### Note!

- *Metal components can be recycled*
- *Upholstery and plastic components must be processed separately*
- *Electrical components (if present) must be disposed of as electronic waste*



#### Warnings

*The product must not be disposed of with household waste!*

### 11.2 Reuse and circularity

To extend the product's service life and support circular use, the product may be offered for reuse or refurbishment.

**To assess the product's suitability for refurbishment or reuse, contact:**

Fysio-Tech B.V. (part of Wesseling B.V.)

## 12. Warranty

### 12.1 Warranty terms

The manufacturer provides a warranty period of 36 months from the date of delivery.

**The warranty covers:**

- Manufacturing defects
- Material defects under normal use

### 12.2 Exclusions

The warranty does not apply to:

- Normal wear of components
- Damage caused by incorrect use or overloading
- Damage caused by external factors (e.g. moisture, impacts)
- Modifications or repairs without the manufacturer's approval

#### Note!

*Use outside the applications described in these instructions may void the warranty*

### 12.3 Repairs

Repairs within the warranty period:

- Are carried out by the manufacturer or an authorised party.

Repairs outside the warranty period:

- Are carried out for a fee.

#### Note!

*Unauthorised repairs may result in unsafe situations and void the warranty.*

---

## 12.4 Spare parts

Use only original spare parts supplied by the manufacturer

### Note!

*Using non-original parts may reduce the safety and performance of the product.*